

Spengler Cardiology presents you *Cardiomate*, a new generation of ECG:

Cardiomate EVI

THE ECG TABLET 100% AUTONOMOUS!



3 TO 5 YEARS GUARANTEE OFFERED
SUPPORT OFFER





With its neat design and compact size, EVI is an intuitive and ultra convenient ECG.

Transportable in all departments and at the bedside of each patient. It is ideal for offices and institutions and meets the requirements of the EN 1789 standard for emergency transport.

Cardiomate EVI has the same performance as the other ECGs of the range *Cardiomate* with a sampling frequency of 32000Hz which gives it a detailed recording of the electrical activity of the heart.

2 sizes of screen 7" and 10" allowing at the same time the tracing on 12 leads as well as the interpretation, and not necessarily requiring printing, for a more ecoresponsible approach.

Cardiomate EVI 10" also exists in 18 leads more useful in emergency settings, in case of chest pain to check for example isolated posterior cardiac attacks, right ventricular heart attacks or other. The *Cardiomate EVI* 10" 18 leads offers 5 waveform display: 18 x1, 12x1 + 6x1, 6x3 +1R, 6x2 + 6x1 and 6x2 + 6x1 + 2R

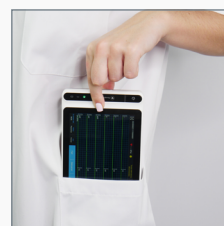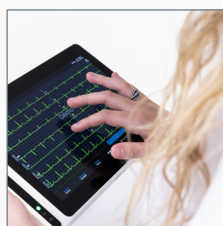

CE
0123
Class IIa



Spengler S.A.S | 30 rue Jean de Guiramand | 13290 Aix en Provence | FRANCE

www.spengler.fr

Tel : +33 (0)4 42 90 31 31
Mail : export@spengler.fr

ECG TABLET

Cardiomate EVI



A TABLET 100% AUTONOMOUS

- 2 sizes : 7" or 10".
- 7" Tablet: ECG 12 leads.
- 10" Tablet: ECG 12 or 18 leads + rechargeable & recorder base.



A FORMIDABLE ACCURACY

- One of the highest sampling rates on the market: 32000Hz.
- Automatic detection and marking of the pacemaker.
- 12-lead ECG or 18-lead ECG possible.



ASSISTANCE SPENGLER

- Pour tous les produits de la gamme EVI, vous bénéficiez d'une assistance gratuite à l'utilisation et au paramétrage de votre tablette ECG Cardiomate EVI. A distance, nos experts vous guident pas à pas.



INTUITIVE AND EFFICIENT

- High resolution 7" or 10" touch screen with keyboard for editing comments on the electrocardiogram.
- Real-time visualization of the 12 ECG leads on the tablet along with automatic interpretation and interval calculations (without using a computer or printer), for a more responsible approach.
- Digital touch signature.



ACCESSIBLE DATA

- Internal memory of 10000 files allowing to store, recall, edit, print patient data.
- Data export via wifi or USB port in pdf or xml format.
- Integration with clinical database systems and hospital management systems.



YOUR FREE DE BOISSY BAG!

- Easily transport your EVI tablet thanks to your De Boissy bag, this historic brand reinvents the traditional doctor's case with style to create practical bags and comfortable that follow you everywhere.



EASY TO PRINT

- Impression of the 12 leads:
- simply on your network printer via wifi for the 2 tablets,
 - or on the thermal printer of the base for the 10".



STAND ALONE

- An ultra-compact design and transportability displayed with:
- protective shell for the 7"
 - desk support for the 10"
 - Compliant with the EN 1789 standard for integration adapted to medical transport vehicles.



A GREAT AUTONOMY

- Both tablets are equipped with a battery of 5000 mAh for an autonomy of at least 12 hours.



www.spengler.fr

Spengler S.A.S | 30 rue Jean de Guiramand | 13290 Aix en Provence | FRANCE

Tel. : **+33 (0)4 42 90 31 31**
Mail : **export@spengler.fr**

ECG TABLET

Cardiomate EVI

TECHNICAL FEATURES

Heart rate range	30bpm~300bpm	Input impedance	≥100MΩ (10Hz)
Accuracy	±1	Input circuit current	≤10nA
Lead	Acquisition of 12-lead Cardiomate EVI 7" and 10" (12 leads) Exists in acquisition of 18-lead Cardiomate EVI 10" (18 leads).	Calibration voltage	1mV±2%
A/D conversion	24 bits	Depolarization voltage	±900mV, ±5%
Sampling rate	32000 samples/Sec	Noise	≤12.5μV
Common mode rejection ratio	≥140dB (AC filter on) ≥120dB (AC filter off)	Amplitude quantisation	0.95μV/LSB
Time constant	≥5s	Recovery time after defibrillation discharge	<10s
Frequency response	0.01HZ~350HZ (+0.4db-3.0db)	Pacer pulse display	Pace beat with an amplitude of ±2mV~±700mV, duration of 0.1ms~2.0ms, A-5 rising time under 100μs, and frequency of 100/min can be displayed on this ECG recording.
Sensitivity	Auto, 2.5 mm/mV, 5 mm/mV, 10 mm/mV, 20 mm/mV, 40 mm/mV, less than ±5% error	Minimum detectable signal	20μVp-p
Filter	AC filter: 50Hz, 60Hz, Off EMG filter: 25Hz, 35Hz, 45Hz, Off ADS filter: 0.01Hz, 0.05Hz, 0.32Hz, 0.67Hz Low pass filter: 75Hz, 100Hz, 150Hz, 300Hz, Off	10.1" Tablet	Size: 272.6mm(L) x 181.2mm(W) x 25.8mm(H) Net weight: 1.2kg
Paper printing Speed	5mm/s, 6.25mm/s, 10mm/s, 12.5mm/s, 25mm/s et 50mm/s, less than ±3% error	10.1" Tablet base	Size: 274mm(L) x 288mm(W) x 111.2mm(H) Net weight: 2.5kg
Analysis algorithm	Glasgow algorithm	7" Tablet	Size: 197mm(L) x 112.4mm(W) x 26.1mm(H) Net weight: 0.8kg
Included	7" Tablet: 1 tablet + shell, + 1 mini medical bag Deboissy Ouessant, 1 patient cable, 4 clamps (terminal electrode), 6 suction balls (chest electrode), 1 power cable, 1 USB key (USB-C USB-A), 1 user manual. 10" Tablet: 1 tablet + charging base with printer, 1 medical bag De Boissy Ouessant, 1 patient cable, 4 clamps (terminal electrode), 6 suction balls (chest electrode), 1 power cable, 1 roll of paper, 1 user manual.		
SPENGLER WARRANTY	3 years including 2 years offered for the 7" tablet 5 years including 4 years offered for the 10" tablet		
COMPLIANCE	EVI tablets comply with the requirements of standard EN 1789 - Medical transport vehicles and their equipment. Road ambulances.		

MODEL	REFERENCE
Cardiomate EVI 7"	260 050
Cardiomate EVI 10" with printing base	260 051
Cardiomate EVI 10" without printing base	260 052
Cardiomate EVI 10" 18 leads (with printing base)	260 054

MODEL	REFERENCE
Medical bag Deboissy Ouessant	MDOUES01
Mini medical bag De Boissy Ouessant	MDOUES02

Accessories

MODEL	REFERENCE
Patient Cable CARDIOMATE EVI 7 & 10"	260 060
Adult limb clamp electrodes for ECG - set of 4	260 070
Pediatric limb clamp electrodes for ECG - set of 4	260 071
Adult suction chest electrode for ECG - set of 6	260 072
Pediatric suction chest electrode for ECG - set of 6	260 073
Roll of paper for CARDIOMATE EVI 10" (20m x 216mm)	260 061
Snap Button Adapters (x10) for CARDIOMATE EVI 7 & 10"	260 023
AC adapter for CARDIOMATE EVI 7 & 10"	260 053

SPENGLER quality

CE Marking

EEC European Directive 93/42: Class IIa. The Cardiomate® meets the essential safety, health, hygiene and environmental protection.

After sales service

Spengler SAS provides maintenance and repair of equipment. Go to the Spengler website in the Contact section.

Spengler Warranty

The Cardiomate EVI is guaranteed 3 years under normal conditions of use.
The warranty does not cover wearing parts, shocks and maintenance defects made by the user.



www.spengler.fr

Spengler S.A.S | 30 rue Jean de Guirmand | 13290 Aix en Provence | FRANCE

Tel. : **+33 (0)4 42 90 31 31**
Mail : **export@spengler.fr**